

OMNIFinger™

Articulating Endoscopic Clip Applier

60° Articulation - 360° Shaft Rotation



Superior Access & Visibility



One-Handed Precision Control



GRENA[®] DG
Think Medical

 **ideal**
Medical Solutions
■ INNOVATIVE EXCELLENCE ■

What is OMNIFinger™?

OMNIFinger™ is an advanced articulating endoscopic clip applier designed to optimize performance in complex laparoscopic procedures. Developed in collaboration with Dr. Diego Gonzalez Rivas, it provides superior access and visibility with 60° of articulation and full 360° shaft rotation.

Purpose-built for demanding surgical scenarios, OMNIFinger™ allows accurate clip placement in hard-to-reach anatomical areas while supporting smooth, ergonomic one-handed operation.

OMNIFinger™ Features

- ✓ **Advanced Articulation:** 60° total angulation with 360° shaft rotation
- ✓ **One-Hand Control:** Index-finger articulation for intuitive operation
- ✓ **Precise Placement:** Enhanced visibility for accurate clip deployment
- ✓ **Versatile Options:** Available in 5mm and 10mm shaft diameters
- ✓ **Flexible Lengths:** 36cm and 47cm working lengths available
- ✓ **Easy Maintenance:** Integrated flushing port for simplified cleaning and reduced contamination risk



Superior Access & Visibility



One-Handed Precision Control




Dr. Diego Gonzalez Rivas
Co-Developer of OMNIFinger™

Dr. Diego Gonzalez Rivas received his Degree in Medicine and Surgery from the Universidad de Santiago de Compostela and completed his residency in Thoracic Surgery at Juan Canalejo Hospital, La Coruña. He has practiced as a Thoracic Surgeon since 2004, with appointments at major university hospitals in Santiago de Compostela and A Coruña.

How OMNIFinger™ Works

OMNIFinger™ is inserted through a compatible 5mm or 10mm trocar, depending on the selected model. With intuitive index-finger control, the distal tip articulates up to 60° while the shaft rotates 360°, enabling precise alignment with the target vessel or structure. Once optimally positioned, compatible Click'aV™ clips are deployed with accuracy and control. The advanced articulating design ensures secure clip placement, even within confined or challenging operative fields.

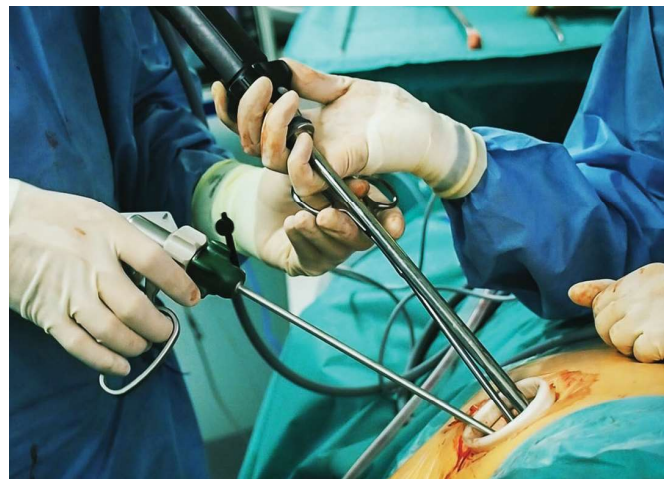

Take a Closer Look

Scan the QR code to explore OMNIFinger™ in more detail and watch the product video.



Ordering Information

REF	SIZE	LOCK	COLOUR	SHAFT DIAMETER (MM)	TROCAR DIAMETER (MM)	WORKING LENGTH (CM)	BOX QUANTITY (PCS)
0301-04MEOMN	M	YES		5	5	36	1
0301-04MEOMNX	M	NO		5	5	36	1
0301-04MLEOMN	ML	YES		5	5	36	1
0301-04MLEOMNX	ML	NO		5	5	36	1
0301-04LEOMN	L	YES		10	10	36	1
0301-04LEOMNX	L	NO		10	10	36	1
0301-04XLEOMN	XL	YES		10	10	36	1
0301-04XLEOMNX	XL	NO		10	10	36	1
0301-04MEOMNB	M	YES		5	5	47	1
0301-04MEOMNBX	M	NO		5	5	47	1
0301-04MLEOMNB	ML	YES		5	5	47	1
0301-04MLEOMNBX	ML	NO		5	5	47	1
0301-04LEOMNB	L	YES		10	10	47	1
0301-04LEOMNBX	L	NO		10	10	47	1
0301-04XLEOMNB	XL	YES		10	10	47	1
0301-04XLEOMNBX	XL	NO		10	10	47	1